

NEW



Video Assisted Heart Valve Surgery



Obadia **II** Instruments Set

Obadia II Instruments Set

DC81115-00

DC81117-94

Delacroix-Chevalier
Exclusive

NEW



**Patent
Pending**

DC81115-00

DC81117-94

Delacroix-Chevalier
Exclusive



Obadia Atrial Retractor for Flexible Additional Blade

NEW

- ✘ Benefit from a full mitral valve exposure thanks to the combination of this new atrial retractor and the flexible additional blade DC81117-94,
- ✘ Keep full advantages of the former Obadia Atrial Retractor,
- ✘ Provided with 3 blades :

35 x 30 mm

35 x 40 mm


35 x 50 mm



DC81115-00 is delivered without short-use flexible additional blade DC81117-94.

DC81117-94

Delacroix-Chevalier
Exclusive

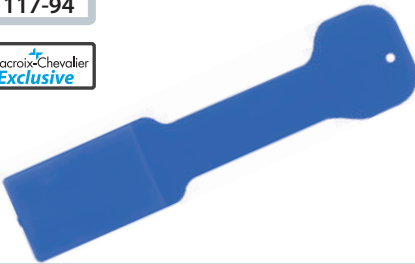

Flexible Additional Atrial Blade

NEW

- ✘ Improves the exposure of the mitral valve by retracting the inferior face of the left atrium and a part of the septum,
- ✘ To be used with the atrial retractor DC81115-00,
- ✘ Short use product (10 uses).

DC81364-36

Delacroix-Chevalier
Exclusive



Forceps for Flexible Additional Atrial Blade

NEW

- ✘ To set flexible additional atrial blade DC81117-94,
- ✘ Can be clipped on ring handle forceps DC81363-24.



DC81363-24

Delacroix-Chevalier
Exclusive



Forceps for Atrial Blade

NEW

- ✘ To hold the atrial blade,
- ✘ Forceps DC81364-36 can be clipped on this forceps for easier gesture.



DC81116-30

DC81116-40

DC81116-50

Delacroix-Chevalier
Exclusive



Spare Blades for Flexible Additional Atrial Blade

NEW

- ✘ Upgrade your Obadia atrial retractor with news blades to benefit from a full mitral valve exposure!
- ✘ To be used with flexible additional atrial blade DC81117-94 and former Obadia atrial retractor DC81110-00.

Obadia II Instruments Set


NEW

Delacroix-Chevalier
Exclusive



DC82152-99

Delacroix-Chevalier
Exclusive

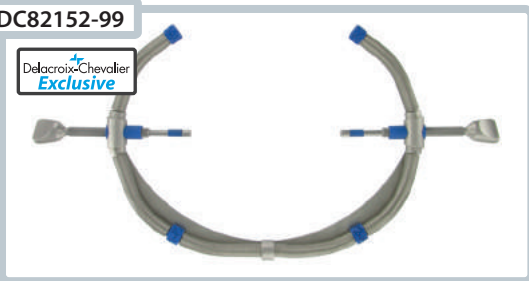

+ Obadia C-Ring Retractor

NEW

- ✘ Get the field free to set the atrial retractor thanks to the opening in the ring,
- ✘ More user friendly due to elliptical shape,
- ✘ For this C-Ring Retractor with stainless steel intercostal blades, ask for DC82152-00:



DC82153-30

Delacroix-Chevalier
Exclusive



+ Soft Intercostal Blades

NEW

- ✘ Protects the tissues of the incision,
- ✘ Retracts the soft tissues,
- ✘ Can provide a light retraction of the ribs thanks to the stainless steel frame,
- ✘ Provided by pair,
- ✘ Short-use product (5 uses).

DC80200-36



MIS DeBakey Forceps - 360 mm - Ø 5 mm

- ✦ Dissecting forceps
- ✦ 2.75mm width jaws



DC80202-36

Delacroix-Chevalier
Exclusive



MIS Atraumatic Resano Forceps - 360 mm - Ø 5 mm

- The Magic Forceps helps during the procedure from sensible to strong grasping
- ✦ Strong grasping atraumatic jaws design



DC80101-36



MIS Curved Scissors - 360 mm - Ø 5 mm

- ✦ Suited for tissues dissection and cutting
- ✦ Not designed for suture cutting



DC80102-36



MIS Strongly Curved Scissors - 360 mm - Ø 5 mm

- ✦ Same benefits as MIS Curved Scissors
- ✦ Suited for difficult access' areas
- ✦ Not designed for suture cutting



DC80012-36

Delacroix-Chevalier
Exclusive



MIS TC Curved Needle Holder - 360 mm - Ø 5 mm

- ✦ New locking system
- ✦ Tungsten inserts for strong and stable grasping
- ✦ From 6/0 to 2/0 suture



DC80019-36

Delacroix-Chevalier
Exclusive



Very Strong MIS TC Curved Needle Holder - 360 mm - Jaws Ø 6 mm

- ✦ New locking system
- ✦ Tungsten inserts for strong and stable grasping
- ✦ From 6/0 to 2/0 suture
- ✦ Very strong tightening



Obadia II Instruments Set

DC81308-35

Delacroix-Chevalier
Exclusive




MIS Knot Pusher - 350 mm

- ✦ Perfect suture guide and retention
- ✦ Push until the very last and critical millimeter
- ✦ A reusable solution



DC81410-28

Delacroix-Chevalier
Exclusive

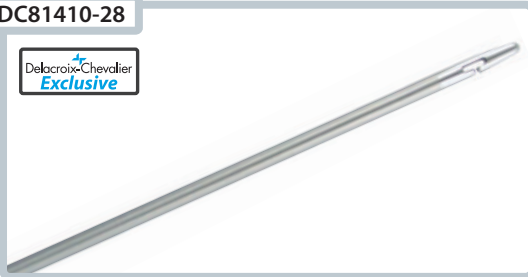

MIS Suture Catcher - 280 mm

- ✦ An exclusive distal tip enabling suture catching through thoracic wall
- ✦ Use for sutures which allow the pericardium exposure



DC81310-28



MIS Magnetic Retriever - 280 mm

- ✦ Could be very useful!
- ✦ An extremely low profile distal tip (Ø: 8,5 mm - L: 8,7 mm)
- ✦ A magnetic power adapted to its use



DC81330-33



MIS Scalpel Blade Holder - 330 mm

- ✦ Easy to set up
- ✦ To use with following blades

Blade Model	Reference
10	B130100
11P	B13011P
12	B130120
15S	B13015S
18	B130180
19	B130190



DC81302-30




MIS Hook for Nerves and Vessel - 300 mm

- ✦ Delicate and smooth surface
- ✦ Precisely laser graduated on request
- ✦ Sabot shape to atraumatically lift and test the Chordae Tendinae




DC81303-33

Delacroix-Chevalier
Exclusive

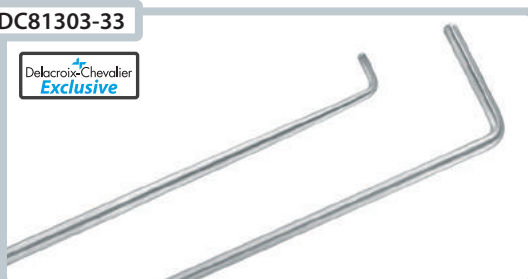

MIS Double Hook for Nerves and Vessels - 330 mm


- ✦ One single instrument for two functions
- ✦ An exclusive sliding handle enabling easy change of configuration and extensive cleaning before sterilization
- ✦ An easy handle locking system
- ✦ One long tip for mitral valve testing

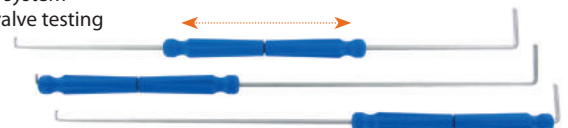

NEW

DC81322-38

Delacroix-Chevalier
Exclusive

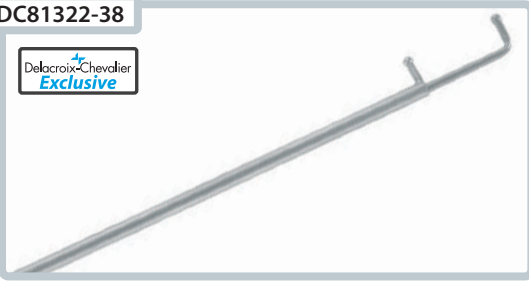

MIS Curved Suture Ruler - 380 mm

- ✦ Essential for chordal replacement procedure
- ✦ Enables to choose the right lengths of artificial chords
- ✦ A knotting device included at distal tip
- ✦ Comes apart for a better handling
- ✦ Exclusive locking system



B25925



Double Curved Scissors - 240 mm

- ✦ User-friendly thanks to double curvature,
- ✦ Strong scissors for the most coarse tissues.




DC80402-32



MIS DeBakey Intercostal Aortic Clamp 2x3 Teeth - 320 mm


- ✦ Its unique shape allows comfortable insertion through 5th inter-costal space, with the efficiency of an endoscopic clamp
- ✦ Proximal part of the clamp pushes left lung and participates to the exposure quality
- ✦ A reusable solution!

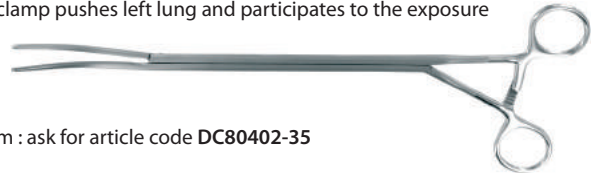

Also available in 350 mm : ask for article code DC80402-35

DC71337-33

Delacroix-Chevalier
Exclusive

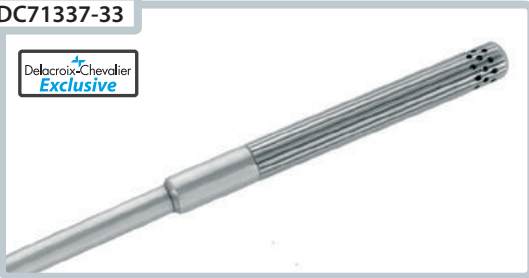

MIS Suction Tube - 330 mm

- ✦ Length adapted to MIS procedures



DC80420-24



Tube Passer - 240 mm

- ✦ Strong grasping jaws
- ✦ To set thoracic drain tubes



Obadia II *Optional Instruments*

DC80207-36



MIS Curved DeBakey Dissecting Forceps - 360 mm

- ✦ Curved jaws for difficult access areas
- ✦ 2.75mm width jaws



DC80110-36

Delacroix-Chevalier
Exclusive



Dedicated Handle for Disposable Scissors Tip - 360 mm Length

- ✦ To be assembled with disposable scissors tips E251105



E251105


Delacroix-Chevalier
Exclusive



Disposable Curved Scissors Tips - Jaws Ø 5 mm

- ✦ Easy assembling with the dedicated handle DC83110-36
- ✦ Perfect cut at every use

DC81320-38

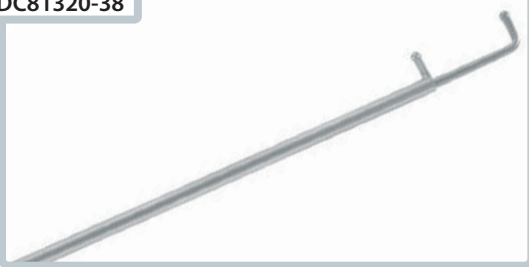


MIS Straight Suture Ruler - 380 mm

- ✦ Essential for chordal replacement procedure
- ✦ Enables to choose the right lengths of artificial chords
- ✦ A knotting device included at distal tip
- ✦ Comes apart for a better handling
- ✦ Exclusive locking system



Obadia II Video Pack Option

E011030WA

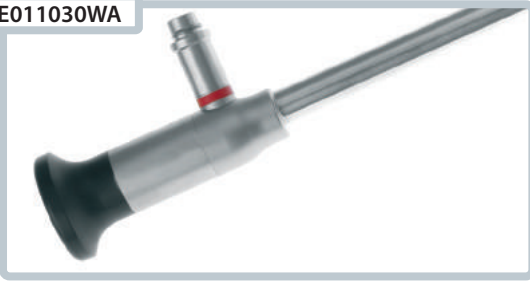

Endoscope Ø 10 mm 30° Angle, HD

- ✦ Superior optic quality, optimized visualization of any cavity
- ✦ Provided with adapter compatible with ACMI/WOLF/STORZ/OLYMPUS



E320001SS-EO




Arm-Strong Endoscope Holder

- ✦ Recommended by Pr. Obadia for fixing the endoscope.
- ✦ Quick and precise position with one hand.
- ✦ Rigid structure for better image quality.
- ✦ Comes with 3 ball joints: Ø 5 mm, Ø 5,5 mm, Ø 10 mm.

Also available in reinforced version. Ask the reference: **E320001SS / S**.

Obadia II Instruments Set

DC81910-01

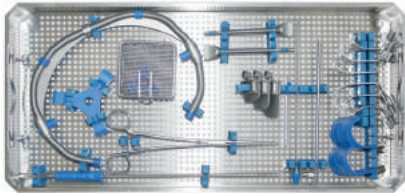

Tray for Obadia II Retractors - 540x255x70 mm

- ✦ For atrial retractor, C-ring retractor and associated blades,
- ✦ Ensures the safety of your instruments during transport, cleaning and sterilization,
- ✦ Recommended to preserve the high level performances of Pr Obadia's instrumentation kit,
- ✦ Provided with lid and silicone brackets.


DC81910-02



Tray for Obadia II Instruments Set - 540x255x70 mm

- ✦ For V handle and round handle instruments,
- ✦ Ensures the safety of your instruments during transport, cleaning and sterilization,
- ✦ Recommended to preserve the high level performances of Pr Obadia's instrumentation kit,
- ✦ Provided with lid and silicone brackets.

DC81910-03

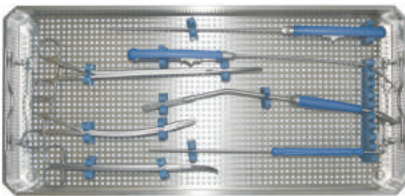

Tray for Obadia II Instruments Set - 540x255x70 mm

- ✦ For aortic clamp, strong scissors, suture rulers and miscellaneous instruments,
- ✦ Ensures the safety of your instruments during transport, cleaning and sterilization,
- ✦ Recommended to preserve the high level performances of Pr Obadia's instrumentation kit,
- ✦ Provided with lid and silicone brackets.

Pr. Jean-François Obadia II Recommended Set

DC81115-00: 3D Atrial Retractor with articulated blade



Included:

DC81116-30: Blade 35 x 30 mm

DC81116-40: Blade 35 x 40 mm

DC81116-50: Blade 35 x 50 mm



DC81117-94: Flexible Additional Atrial Blade



DC81364-36: Forceps for Flexible Additional Atrial Blade



DC81363-24: Forceps for Atrial Blade



DC80200-36: MIS DeBakey Forceps - 360 mm



DC80202-36: MIS Atraumatic Resano Forceps - 360 mm



DC80101-36: MIS Curved Scissors - 360 mm



DC80102-36: MIS Very Curved Scissors - 360 mm



DC80012-36: MIS TC Curved Needle Holder - 360 mm



DC80019-36: Very Strong MIS TC Curved Needle Holder - 360 mm



DC81308-35: MIS Knot Pusher



DC81410-28: MIS Suture Catcher



DC82152-99: C-Ring Retractor



DC82153-30: Soft Intercostal Blade (pair)



DC81310-28: MIS Magnetic Retriever



DC81330-33: MIS Scalpel Blade Holder



DC81302-30: MIS Hook for Nerves and Vessel



DC81303-33: MIS Double Hook for Nerves and Vessel



DC81322-38: MIS Curved Suture Ruler - 380 mm



B25925: Double Curved Scissors - 240 mm



DC80402-32: MIS DeBakey Intercostal Aortic Clamp - 320 mm



DC71337-33: MIS Suction Tube - 330 mm



DC80420-24: Tube Passer - 240 mm



DC81910-01: Sterilization Trays - 540x255x70 mm

DC81910-02: Sterilization Trays - 540x255x70 mm

DC81910-03: Sterilization Trays - 540x255x70 mm



Pr. Jean-François Obadia II Optional Instruments

DC80207-36: MIS Curved DeBakey Dissecting Forceps



DC80110-36: Dedicated Handle for Disposable Scissors Tips



DC81320-38: MIS Suture Ruler



E251105: Disposable Curved Scissors Tips



Pr. Jean-François Obadia II Video Pack

E011030WA: Endoscope 30° Angle HD



E320001SS-EO: Arm-Strong Endoscope Holder




Delacroix-Chevalier
Paris

Corporate Office

13. Avenue de la République - 75011 Paris - France

Manufacturing and logistic

Rue Renée et Louis Landanger

52000 Chaumont - France

Phone : +33 (0) 3 25 02 10 35 - Fax : +33 (0) 3 25 02 10 20

Email : export@delacroix-chevalier.com www.delacroix-chevalier.com

© Copyright 2009-2017 / The photos and pictures are not contractual. Patterns are subject to change. / ZA-B00002